

Telephone dialler ProfiCaller

The **ProfiCaller** telephone dialler independently calls up to 4 telephone numbers in the event of a fault and issues a previously entered text. There is a free choice of dialling sequence. The automatic redial setting is 3 times for each number. The message text has a maximum length of 40 seconds. The dialler comes with a plug-in power supply unit. All 4 alarm outputs are preconfigured in the standard (5 m connecting lead).



ProfiCaller

External Dimensions (WxDxH) mm	: 150x30x104
Mains connection volts/Hz/ph	: 230/50/1 with a plug-in power supply unit 11.5V up to 14V
PBX-capable	: yes
System requirements	: telephone connection (TAE-N socket); no-volt remote alarm contact at the cooling and freezing devices or Temperature Monitoring and Recording System ProfiTec (see separate brochure)

ProfiCaller + Accessories:

Cat. No.

ProfiCaller	ZXX2-009-TX
Uninterrupted power supply unit (UPS).	ZXX2-007-20


Connection options:

For up to 4 cooling and freezing devices **with** no-volt remote alarm contact.

- ① Telephonic alarm via ProfiCaller.


For up to 4 cooling and freezing devices **without** no-volt remote alarm contact.

- ② Acoustic and optical alarm via ProfiTec (see separate brochure).
- ③ Telephonic alarm via ProfiCaller via ProfiTec.

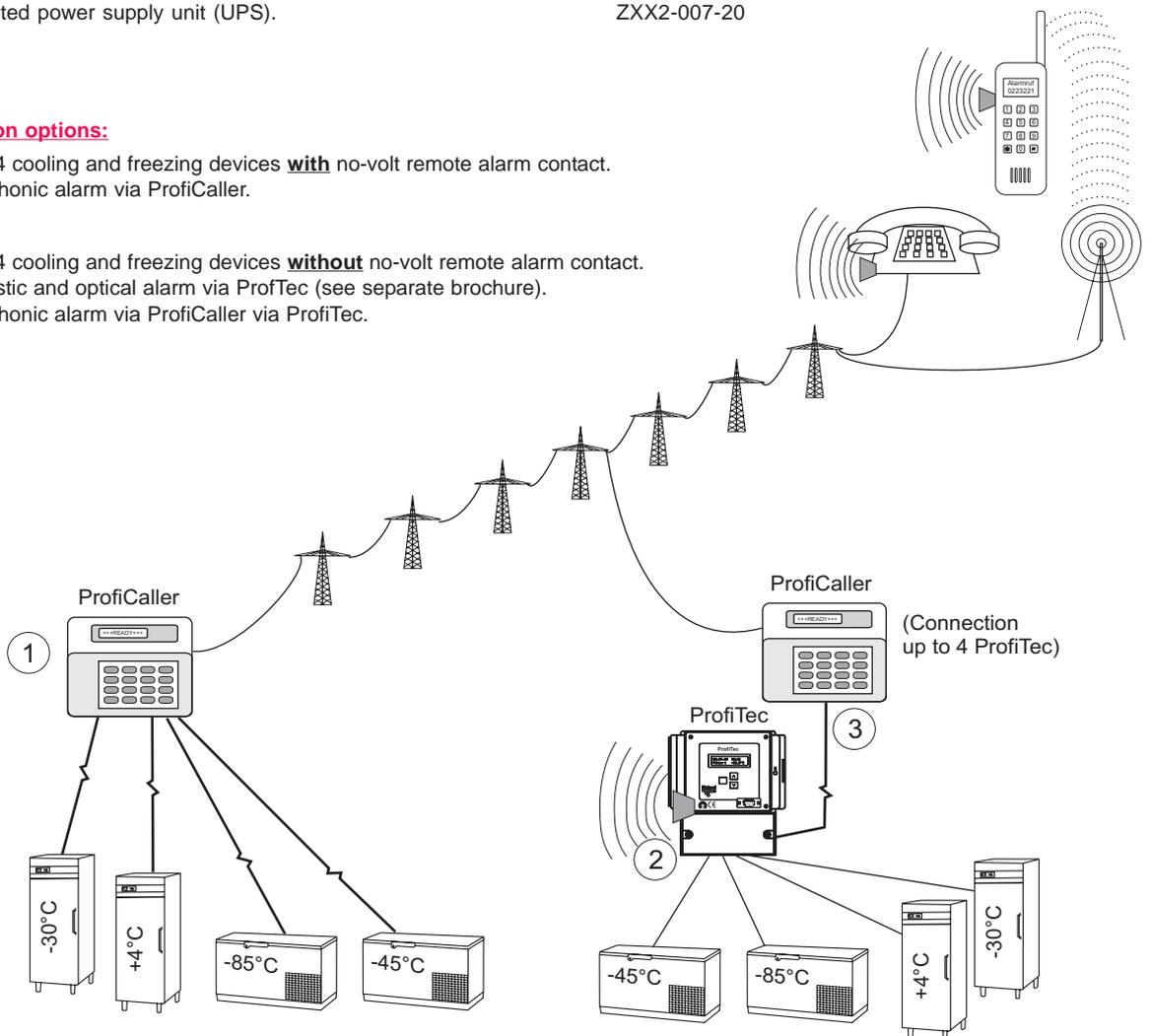

Chilling or freezing device **with** no-volt remote alarm contact

Chilling or freezing device **without** no-volt remote alarm contact

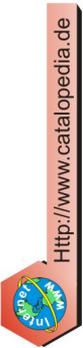

Http://www.catalopedia.de
 e-Mail: LabSales@NationalLab.com

National Lab GmbH
 Grambeker Weg 157
 D-23879 Mölln

Tel.: +49 45 42 / 84 91 - 70
 Fax: +49 45 42 / 84 91 - 71